# BabyKit

**EN – MANUAL AND SAFETY INSTRUCTIONS**

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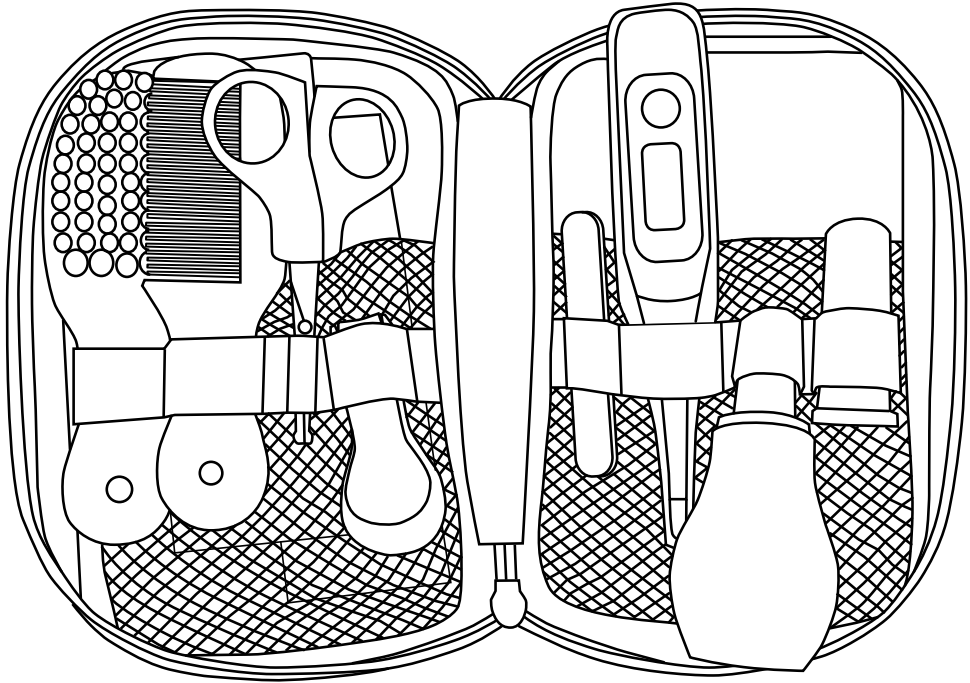
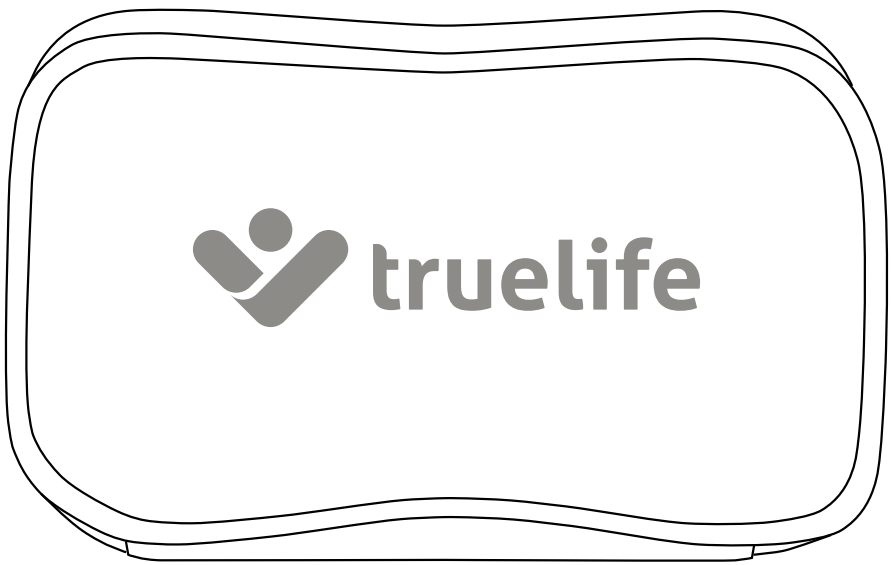
---

### PACKAGE CONTENTS

- |                                  |  |
|----------------------------------|--|
| 1× mucus extractor               | 3× nail file   |
| 1× round nail scissors           | 1× manual and safety instructions  |
| 1× nail clippers                 | 1× TrueLife Care T3 thermometer (instructions and safety instructions for the thermometer separately in the package) |
| 1× comb                          |  |
| 1× hair brush with fine bristles |  |
| 1× finger toothbrush             |  |


### MUCUS EXTRACTOR

#### Notice:

- Do not push the mucus extractor deep into the nasal cavity, there is a risk of injury.
- Under no circumstances should you insert the mucus extractor into body openings other than the nose.
- Do not use the extractor more than 4 times a day, there is a risk of irritation of the child's mucous membrane.
- Use the extractor slowly and gently, rapid and repeated use may cause irritation or inflammation of the nasal tissue.

#### Procedure for use:

- 1) Squeeze the mucus extractor and carefully insert it into the child's nostril.
- 2) Slowly release the pressure on the suction cup and remove it from the nose.
- 3) Open the aspirator, empty the mucus into a tissue and close it again.
- 4) Repeat the same procedure for the other nostril.
- 5) Disassemble the extractor, rinse it in hot water and then dry it.
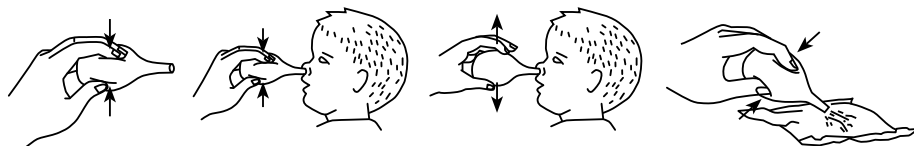
#### Cleaning:

- Disinfect the mucus extractor before each use to prevent the transfer of bacteria to the child.
- Suitable disinfection methods are:
  1. Immersion in boiling water for at least 15 minutes.
  2. Use of a disinfectant gel designed for infant needs.

### FINGER TOOTHBRUSH

- 1) Put the brush on your finger.
- 2) Carefully, with gentle movements, run over the child's gums and along the tongue. For the smallest children, cleaning the oral cavity should not take more than 30 seconds.
- 3) Rinse the brush with hot water and then let it dry.

#### Cleaning:

- Disinfect the toothbrush before each use to prevent the transfer of bacteria to the child.
- Appropriate disinfection methods are:
  1. Immersion in boiling water for at least 15 minutes.
  2. Use of a disinfectant gel designed for infant needs.

### ROUND NAIL SCISSORS AND NAIL NIPPERS

#### Notice:

- Never cut the nails too close to the baby's skin. Risk of injury!

- We recommend trimming the child's nails after a bath, when they are softened by water and easier to cut, or when the child falls asleep.

**Procedure for use:**

- 1) Grasp the baby's finger between your index finger and thumb and use the rest of your hand to cover the baby's other fingers.
- 2) Carefully trim the nail.
- 3) Repeat the procedure for other fingers.

**Cleaning:**

- You can clean the blades of scissors and nail clippers with a wet tissue or disinfectant gel designed for baby supplies.

**NAIL FILES****Procedure for use:**

- 1) Slowly and carefully run the file over individual nails.

**HAIR BRUSH WITH FINE BRISTLES AND COMB****Notice:**

- Avoid pulling the baby's hair. When brushing or combing, hold the roots of the hair with your hand so as not to pull the child's scalp unnecessarily.

**Procedure for use:**

- 1) Brush and comb baby's hair slowly and carefully. Work from roots to ends.

Before first use, the consumer is obliged to familiarize himself with the principles of safe use of the product.

**POLICIES AND NOTICES**

- Please read these safety instructions first before using the product.
- For proper handling of the product, read the most current version of the manual, which can be downloaded from the product page at: <https://www.truelife.eu>.

**MORE INFORMATION**

Apply for device warranty repairs at your dealer. In case of technical problems and questions, contact your dealer, who will inform you about the next procedure.

The warranty period for products is 24 months, unless otherwise specified. The warranty does not cover damage caused by non-standard use, mechanical damage, exposure to aggressive conditions, handling contrary to the manual and normal wear and tear. The warranty period for the battery is 24 months, for its capacity 6 months. More information about the warranty can be found at [www.elem6.com/warranty](http://www.elem6.com/warranty)

The manufacturer, importer or distributor shall not be liable for any damage caused by the installation or misuse of the product.

Typographical errors and changes in the manual and safety instructions are reserved.

**Importer:**

elem6 s.r.o.,  
Braškovská 15, 16 100 Prague 6, Czech Republic  
Made in P.R.C.

JOYTECH Healthcare Co., Ltd.  
No.365, Wuzhou Road, Yuhang Economic Development Zone,  
Hangzhou City, 311100 Zhejiang P.R. China

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